

# Herbimycin A Inhibitor

## Catalog# SIH-116A/B

Size: 100ug/1mg

PO Box 30244, Suite 405,  
3989 Quadra Street,  
Victoria, BC V8X 5E1, Canada

This product is for *in vitro* research use only and is not intended for use in humans or animals

# StressMarq

## Biosciences Inc.

Orders • [sales@stressmarq.com](mailto:sales@stressmarq.com)

Tel: • +1 250 294 9065

Fax: • +1 250 294 9025

Email • [info@stressmarq.com](mailto:info@stressmarq.com)

Web • [www.stressmarq.com](http://www.stressmarq.com)

Product	Herbimycin A
Formula	C <sub>30</sub> H <sub>42</sub> N <sub>2</sub> O <sub>9</sub>
MW	574.3
Source/Host	Synthetic
Purity	98%
Solubility	Soluble in DMSO and ethanol (1mg/mL)
Appearance	Lyophilized solid; yellow
Storage and stability	-20°C; 1 year+; shipped ambient Protect from light.

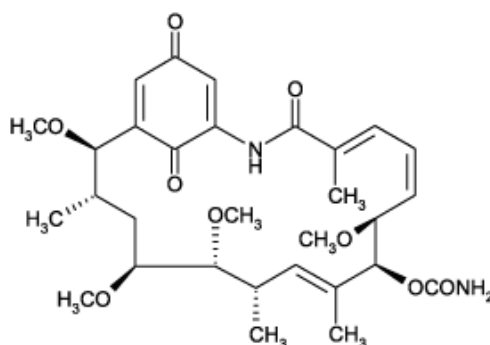


Figure 1: Structure of Herbimycin

### Scientific Background

Herbimycin A, a benzoquinoid ansamycin antibiotic, irreversibly and selectively inhibits tyrosine kinases by reacting with thiol groups (1,2). It is effective on Src, Yes, Fps, Ros, Abl and ErbB oncogene products (1-3) and inhibits thrombin-induced tyrosine phosphorylation of phospholipase C (4).

### Selected References

1. Uehara Y., and Fukazawa H. (1991) *Methods Enzymol.* 201: 370.
2. Fukazawa H., *et al.* (1991) *Biochem. Pharmacol.* 42: 1661.
3. Satoh T., *et al.* (1992) *J. Biol. Chem.* 267: 2537.
4. Weiss R., and Nuccitelli R. (1992) *J. Biol. Chem.* 267:5608.

# Material Safety Data Sheet

This product is for *in vitro* research use only and is not intended for use in humans or animals

The below information is believed to be correct but does not purport to be all inclusive and shall be used only as a guide. StressMarq shall not be held liable for any damage resulting from handling or from contact with the above product. See the Technical Specification, Packing Slip, Invoice, and Product Catalogue for additional terms and conditions of sale.

---

## Hazardous Ingredients

The physical, chemical and toxicological properties of these components have not been fully investigated. It is recommended that all laboratory personnel follow standard laboratory safety procedures when handling this product. Safety procedures should include wearing OSHA approved safety glasses, gloves and protective clothing. Direct physical contact with this product should be avoided.

<u>Known Hazardous Components</u>	<u>CAS Number</u>	<u>Percent</u>
None		

---

## Physical Data

This product consists of powder shipped at ambient temperatures. The physical properties of this product have not been investigated thoroughly. CAS number 70563-58-5, Chemical Class: Macrolide

---

## Fire and Explosion Hazard and Reactivity Data

NOT APPLICABLE

---

## Toxicological Properties

May be harmful by inhalation, ingestion, or skin absorption. The toxicological properties of this product have not been investigated thoroughly. Exercise due caution.

---

## Preventative Measures

Wear chemical safety goggles and compatible chemical-resistant gloves. Avoid inhalation, contact with eyes, skin or clothing.

---

## Spill and Leak Procedures

Observe all federal, state and local environmental regulations.

- Wear protective equipment.
- Absorb on sand or vermiculite and place in closed containers for disposal.
- Dispose or mix the material with a combustible solvent and burn in a chemical incinerator equipped with an afterburner and scrubber.

---

## First Aid Measures

- If swallowed, wash out mouth with water, provided person is conscious. Call a physician.
- In case of skin contact, flush with copious amounts of water for at least 15 minutes. Remove contaminated clothing and shoes. If a rash or other irritation develops, call a physician.
- If inhaled, remove to fresh air. If breathing becomes difficult, call a physician.
- In case of eye contact, flush with copious amounts of water for at least 15 minutes while separating the eyelids with fingers. Call a physician.

Authorized: StressMarq Biosciences Inc.  
Creation Date: 06/20/2008